



MIM

METAL INJECTION MOLDING (MIM)

can produce small, highly complex parts with tolerances as small as 0.0005 inches. MIM's advantages lie in its capability to produce mechanical properties with good dimensional tolerance control.

Metal injection molded parts offer a nearly unlimited shape and geometric-feature capability. Lancer's brackets are created to meet exact design specifications for a robust, high-quality, dimensionally accurate appliance.



Milling

Because nothing but the highest levels of accuracy will do, Lancer uses a precision milling process to form exact archwire slots. Milled brackets are dimensionally accurate and highly functional, allowing for the fullest expression of the prescription.




BEFORE




MILLING



AFTER

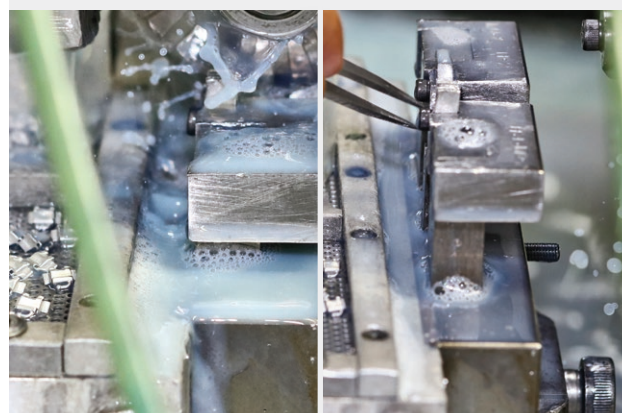

PRAXIS™ DELIVERS THE HIGHEST LEVEL OF PRECISION

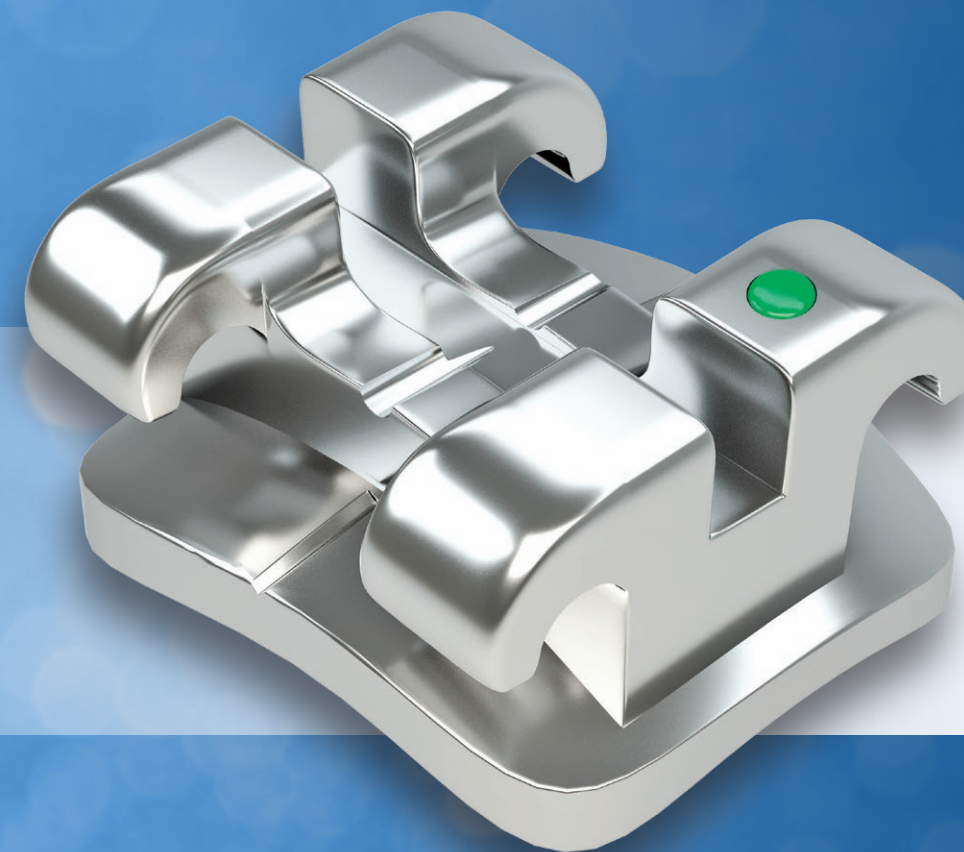

LANCER
ORTHODONTICS



EST. 1967

LANCER
ORTHODONTICS



EST. 1967

(800) 854-2896 • www.lancerortho.com


International: +1 (760) 744-5585 • Fax: +1 (760) 598-0418
2726 Loker Avenue West • Carlsbad, California 92010 (U.S.A.)
info@lancerortho.com

Treat with High Precision

Praxis™ Mini-Twin

A perfect combination of robotic MIM technology, superior design, and pleasing aesthetics, producing a robust bracket that is inherently stronger and more dimensionally accurate than similar brackets.

This powerful, functionally superior bracket is loaded with features that patients and doctors love.

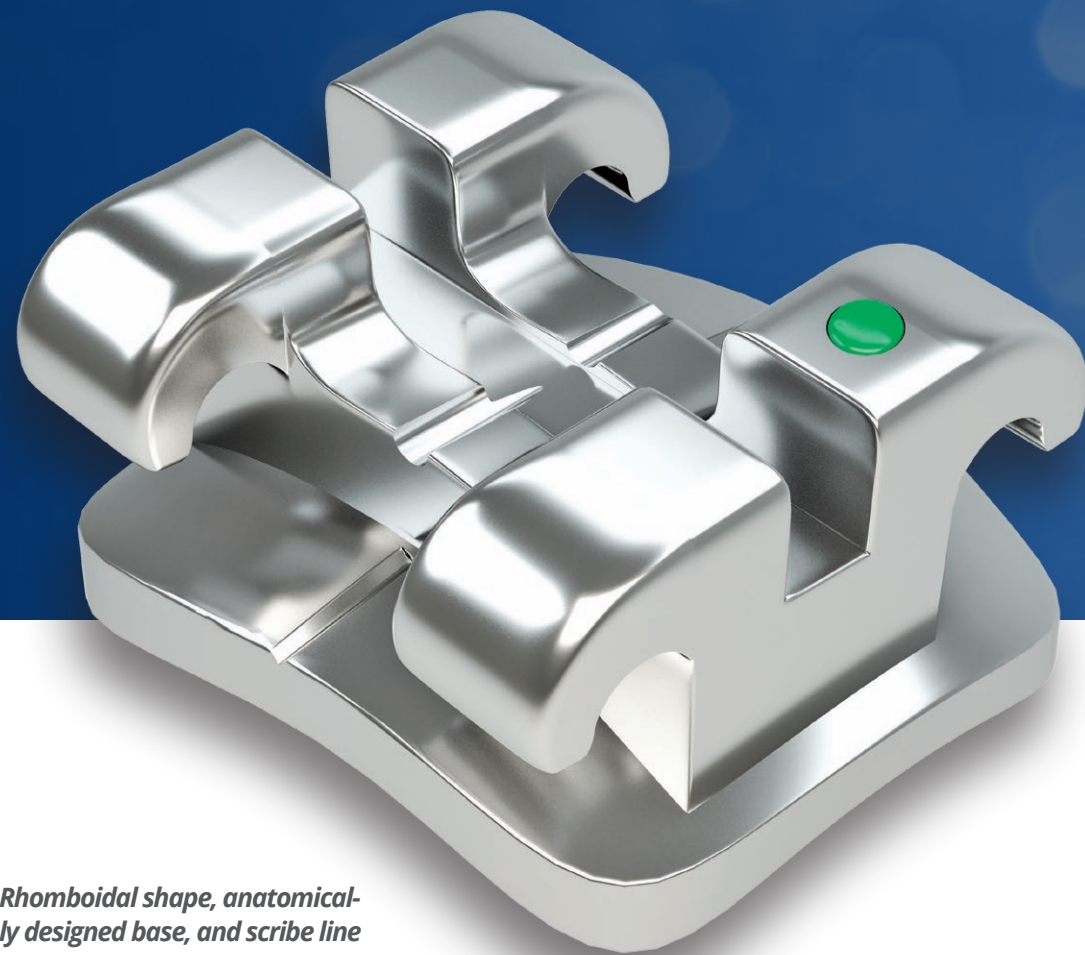

MILLED



MIM



MESH BASE



REDI-PAK™

- Rhomboidal shape, anatomically designed base, and scribe line which runs through both the bracket and the pad facilitate fast and accurate placement

- Slots are milled with diamond blades, guaranteeing precise dimensions throughout the wall of the slot

- Torque-in-base allows for less wire bends and fewer brackets being repositioned

- Low-profile design for optimum patient comfort

- Large under tie wing area for easy ligation

- Anatomically designed 80-gauge mesh pads make bonding effortless while ensuring optimal adhesion

LANCER'S REDI-PAK™ COMES IN A VARIETY OF PRESCRIPTIONS

and keeps brackets in individualized, sealed packages until patients are ready for bonding. When it's time to treat, unfolding the Redi-Pak puts the brackets in "tweezer ready" position. And, if storage is required, Redi-Pak can be resealed and placed with your patient's chart.



Roth Redi-Pak™

		Slot .018	Slot .022
	Hooks		
All Hooks	x	265-933	265-934
Cuspid Hooks	x	265-937	265-938
No Hooks		265-941	265-942



MBT Redi-Pak™

		Slot .018	Slot .022
	Hooks		
All Hooks	x	268-931	268-932
Cuspid Hooks	x	268-933	268-934
No Hooks		268-935	268-936

Also available in Hilgers, Ricketts and Andrews Prescriptions.

Praxis Gold™

Offer a stylish alternative and maximize patient satisfaction with Praxis Gold™.

Precisely engineered with all the classic features of our traditional Praxis™ brackets, Praxis Gold™ offers an elegant 24 karat gold coating that will maintain its brilliant luster throughout the duration of treatment without chipping or peeling.

Safe for use on nickel-sensitive patients, the Praxis Gold™ is a stylish alternative to traditionally finished brackets.



Praxis Gold™ MBT Redi-Pak™

		Slot .022
All Hooks	x	268-992

